

Precipitation and Water Purity

SJTU JI


&

Department of Chemistry, SJTU

Course Evaluation!
5 min

Goals for Experiment

- Part 1: precipitation rxn & filtration
- Part 2: Designated cations groups I & II precipitation trends with designated anions
- Part 3: Effect of cations & anions concentration on precipitation. Q_{sp} trends relative to K_{sp}
- Part 4: Solubility of salts with polar & non-polar solvents: I. CaCl_2 & II. $\text{K}_2\text{C}_2\text{O}_4$ solubilities in H_2O , Acetone, & Hexane.
- Part 5: Omit.
- Complete report during lab work.

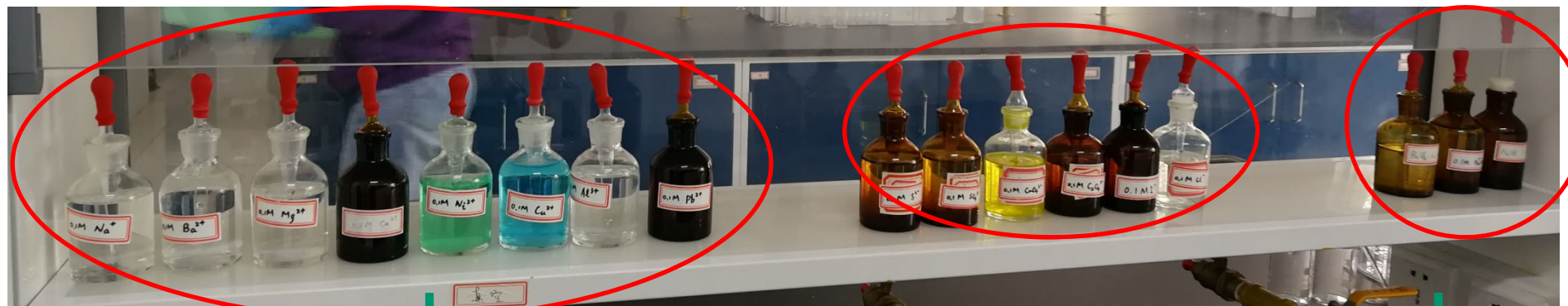


Cations (8)

Anions (6)

BaCl_2
 NaOH

6 Dry tubes & Filtration paper @ Front Desk
Acetone & Hexane @ Side Countertops
 CaCl_2 & $\text{K}_2\text{C}_2\text{O}_4$ @ Side Countertops

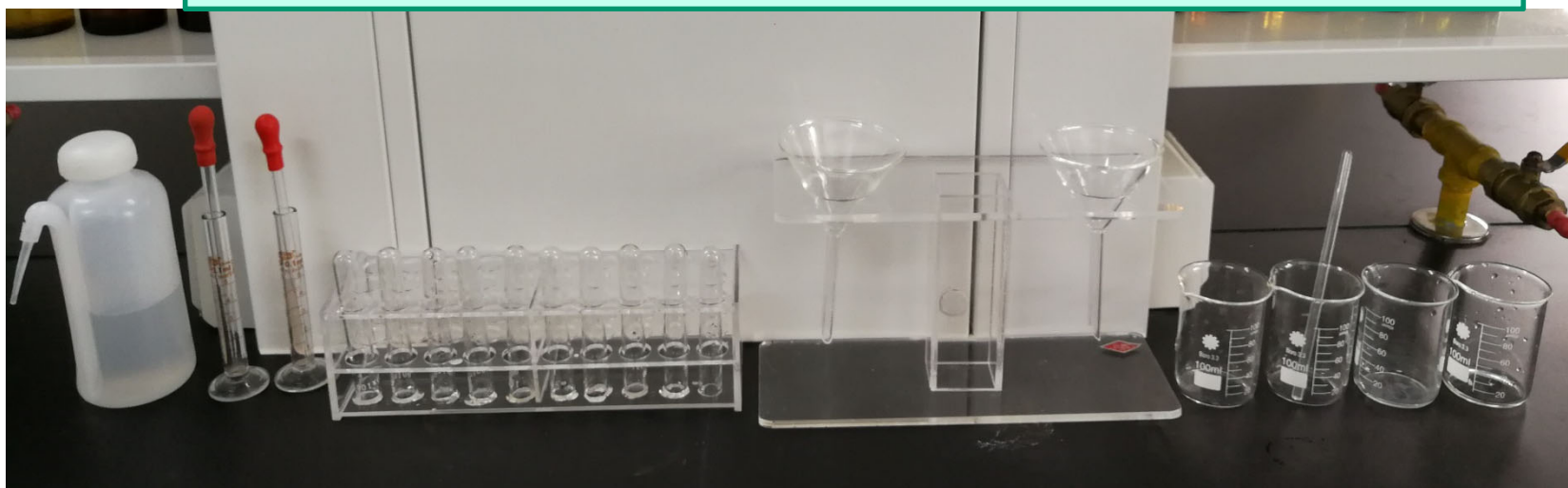

Part 1: What is a Precipitate?

CuSO_4 & BaCl_2
5 mL each in Beaker




Settle >30 min
NO stirring!



Filtration
Gravity

Demo

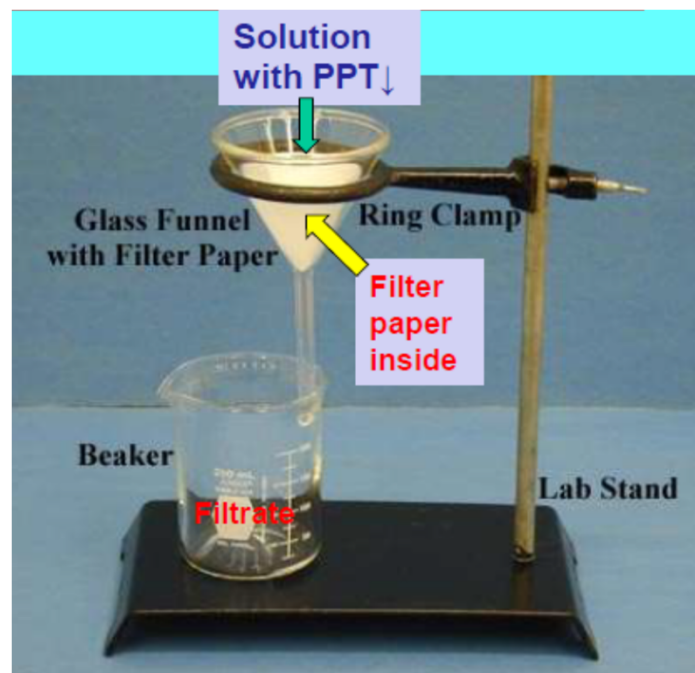

2 Stud.

Part 2: Predictable?

Cation & Anion

2 drops in tubes



Use the A4 Table

Record!!!

Tabl1: CATION GROUPS TO PRECIPITATE (0.10 M nitrate salts for each cation), p52 manual.

GROUP I	Na ⁺	Ba ²⁺	Mg ²⁺	Co ²⁺	Ni ²⁺	Cu ²⁺	Al ³⁺	Pb ²⁺
GROUP II	K ⁺	Mn ²⁺	Ca ²⁺	Sr ²⁺	Cr ³⁺	Fe ³⁺	Zn ²⁺	Ag ⁺

Table 2: ANION GROUPS PRECIPITATING REAGENTS: (0.10 M sodium salts for each anion), p52 manual. SELECT CATION GROUP I OR GROUP II BUT NOT BOTH, THEN ADD YOUR OBSERVATIONS ON THE TABLE ACCORDING TO INSTRUCTIONS & USING GROUP II EXAMPLE for Ag⁺ (CAUTION: COLORS MAY NOT BE AS SHOWN). COMPARE YOUR RESULTS WITH THE OTHER GROUP I/II. Clear means no ppt.

CATION GROUP I 2drops	CATION GROUP II 2drops	REF 2drops	Cl ⁻ 2drops	CrO ₄ ²⁻ 2drops	I ⁻ 2drops	C ₂ O ₄ ²⁻ 2drops	S ²⁻ 2drops	SO ₄ ²⁻ 2drops	SPCTR GROUP I	SPCTR GROUP II
Na ⁺	K ⁺	clear								
Ba ²⁺	Mn ²⁺	clear								
Mg ²⁺	Ca ²⁺	clear								
Co ²⁺	Sr ²⁺	clear								
Ni ²⁺	Cr ³⁺	clear								
Cu ²⁺	Fe ³⁺	clear								
Al ³⁺	Zn ²⁺	clear								
Pb ²⁺	Ag ⁺	clear	White ppt↓	Brown ppt↓	Yellow ppt↓	White ppt↓	Black ppt↓	White ppt↓		

Odd #: Group I
Even #: Group II

4 Stud.

Part 3: Conc. & Precip.

Design

- 3 Cations
 - NaOH & 1 other anion
(1 drop w/ or w/o 9 drops of DI)
- In Tubes



Record!

Reaction #	Reactant #1	Reactant #2
II-1	0.10M Pb(NO ₃) ₂	0.10M NaOH
II-2	0.10M Pb(NO ₃) ₂	1.0M NaOH
II-3	0.01Pb(NO ₃) ₂	0.01M NaOH
IV-1	0.10M ZnSO ₄	0.10M NaOH
IV-2	0.10M ZnSO ₄	1.0M NaOH
IV-3	0.01M ZnSO ₄	0.01M NaOH
V-1	0.10M CaCl ₂	0.10M K ₂ C ₂ O ₄
V-2	0.01M CaCl ₂	0.01M K ₂ C ₂ O ₄

4 Stud.




Part 4: Solvent

Aqueous
Organic (dry tubes)
Solid @ side countertop
~2 mL by counting drops



Record!

- Water (H_2O) is very polar.
- Hexane (C_6H_{14}) is nonpolar.
- Acetone (CH_3COCH_3) is moderately polar.

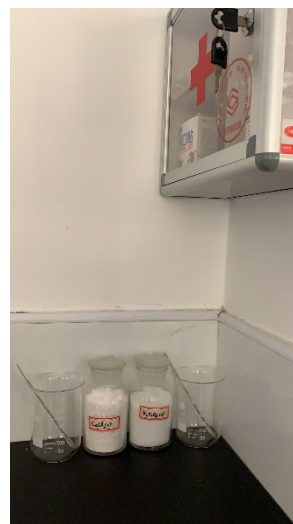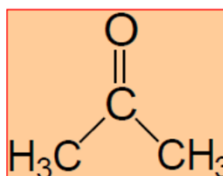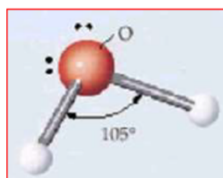

4 Stud.

Report

PLE: Pre-lab Report (p18-19) Individual

- Part 2.A., A4 Table w/ plastic cover
(Hypothesis and predications for question 4 can be included in the A4 table)
- Prepare your own data sheet & team assessment form

PLQ: Team Report (p20-22) Group & Individual

- Omit Part 2.A. Data Analysis (p21)
- Omit Part 2.B. (p21)
- Omit Part 5 (p22)
- 'class data', two groups working on the two sides of one counter complete a whole test set and the data can be considered as class data.
- Data sheet will be graded individually

ALR: Individual Report Individual

- 4 parts: Intro., Procedures, Discussion (p23 1-5), and Reference.
- Procedure, briefly summarize the operation steps.