

**SAFETY WARNING: CAUSTIC CORROSIVE TOXIC CHEMICALS & PRESENCE OF FLAME/HEAT.**

1. **YOU MUST WEAR SAFETY GOGGLES ALL THE TIME.**
2. **MUST HANDLE THE CHEMICALS & HEATING EQUIPMENT VERY CAREFULLY, AVOID CONTACT WITH SKIN, RINSE WITH PLENTY OF WATER IMMEDIATELY UPON CONTACT.**
3. **MUST DISPOSE WASTE CHEMICALS AS DIRECTED & PLACE IN THE APPROVED WASTE CONTAINERS.**

**4. YOU MUST OBEY ALL THE SAFETY RULES OR YOU MUST LEAVE THE LABORATORY IMMEDIATELY.**

**\* DO NOT WASTE CHEMICALS.**

**\* DO NOT REPEAT TRIALS MORE THAN INSTRUCTED.**

**\* DO E4(I) EXPERIMENT FIRST, EACH TWO STUDENTS IN A GROUP MUST DO THE ENTIRE EXPERIMENT AS INSTRUCTED BY E4(I).**

**\* DO E4(II) EXPERIMENT NEXT AS INSTRUCTED BY E4(II). EACH GROUP DO THE ENTIRE EXPERIMENT WHERE EACH INDIVIDUAL IN A GROUP MUST DO 1 SAMPLE REACTION REPEATED ONCE ONLY.**

**\* SHARE RESULTS WITH YOUR GROUP AND COMPLETE INDIVIDUAL REPORTS.**

***Raw Data Sheet!!!***

**Introduction to Kinetics:**  
**Factors That Affect the Rate of Reaction**  
**Determining the Rate Law:**  
**A Kinetics Study of the Iodination of Acetone**

**Objectives:**

1. Be able to list and rationalize the factors that affect the rate of reaction.
2. Explain various scenarios using the factors that affect reaction kinetics.
3. Gain a quantitative understanding of kinetics.
4. Determine the rate of a reaction, the order of the reaction with respect to the reactants and the value of the rate constant.
5. Predict reaction times using an experimentally results of the rate law.

# Chemical Kinetics

*Study of the rates of chemical reactions*

Reaction rates vary greatly

**very fast** (e.g., combustion)

**very slow** (e.g., disintegration of a plastic bottle in sunlight).

***Concentrations of reactants***

***Temperature***

***Surface area***

***Catalysts***

# Determining the Rate Law:

## A Kinetics Study of the Iodination of Acetone

**The Rate Law for RXN:**  $aA + bB \rightarrow cC + dD$

Rate of Disappearance of A = Reaction Rate =

$$-R_A = d[A]/dt = v = k[A]^m[B]^n$$

**Where m and n** are determined experimentally.

**k** = Reaction Rate Constant (its units depends on reaction).

**Arrhenius Equation:**  $k = f e^{-E_a/RT}$

Where **f** is frequency factor, **E<sub>a</sub>** is activation energy, **R** is gas constant, and **T** is the absolute temperature

# Rate Laws

The *values of the exponents* in a rate law establish the **order of a reaction**

For reactant A,

$$\text{Rate} = k[A]^m[B]^n$$

if  $m = 1$ , reaction is first order in A

if  $m = 2$ , reaction is second order in A

- The proportionality constant,  **$k$ , is the rate constant** and its value depends on the reaction, the temperature, and the presence or absence of a catalyst.
- Reaction rate constant exponentially proportional<sup>EOS</sup> to temperature (Arrhenius Equation).

# Introduction to Kinetics: Factors That Affect the Rate of Reaction

## A. Effect of Changing the Concentration of Reactants:

<b>Solid:</b>	<b>Effect of Changing the Surface Area</b>
<b>Gas:</b>	<b>Pressure of the Gas</b>
<b>Liquid:</b>	-----
<b>Solution:</b>	<b>Concentration of the Reactant</b>

## Solid: Effect of Changing the Surface Area

**Gas:** **Pressure of the Gas**

**Liquid:** -----



<b>Solution:</b>	<b>Concentration of the Reactant</b>
------------------	--------------------------------------

# Introduction to Kinetics:


## Factors That Affect the Rate of Reaction

### B. Effect of Changing the Temperature:

(According to thermodynamics for reactions at equilibrium,  $K_{\text{equil}}$  decreases with increased temperature so the rate of reaction from thermodynamics view decreases). **See next slide.**

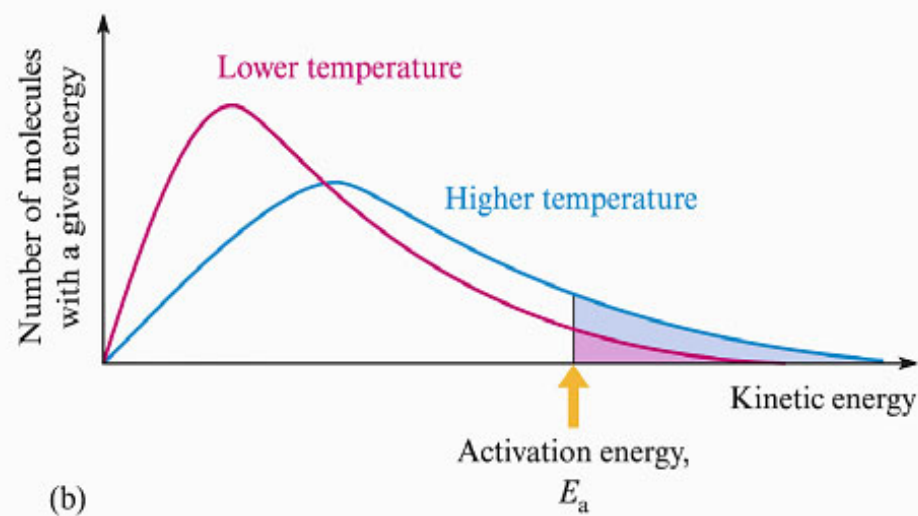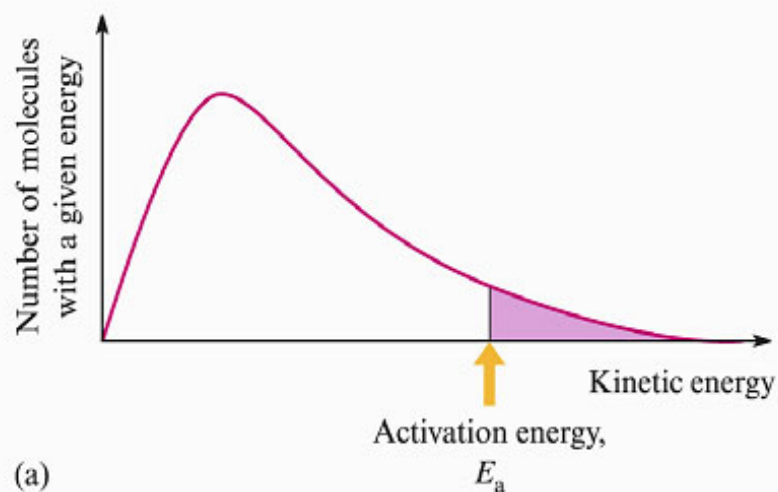

# AUTOMOBILE CATALYTIC CONVERTER OPTIMUM DESIGN

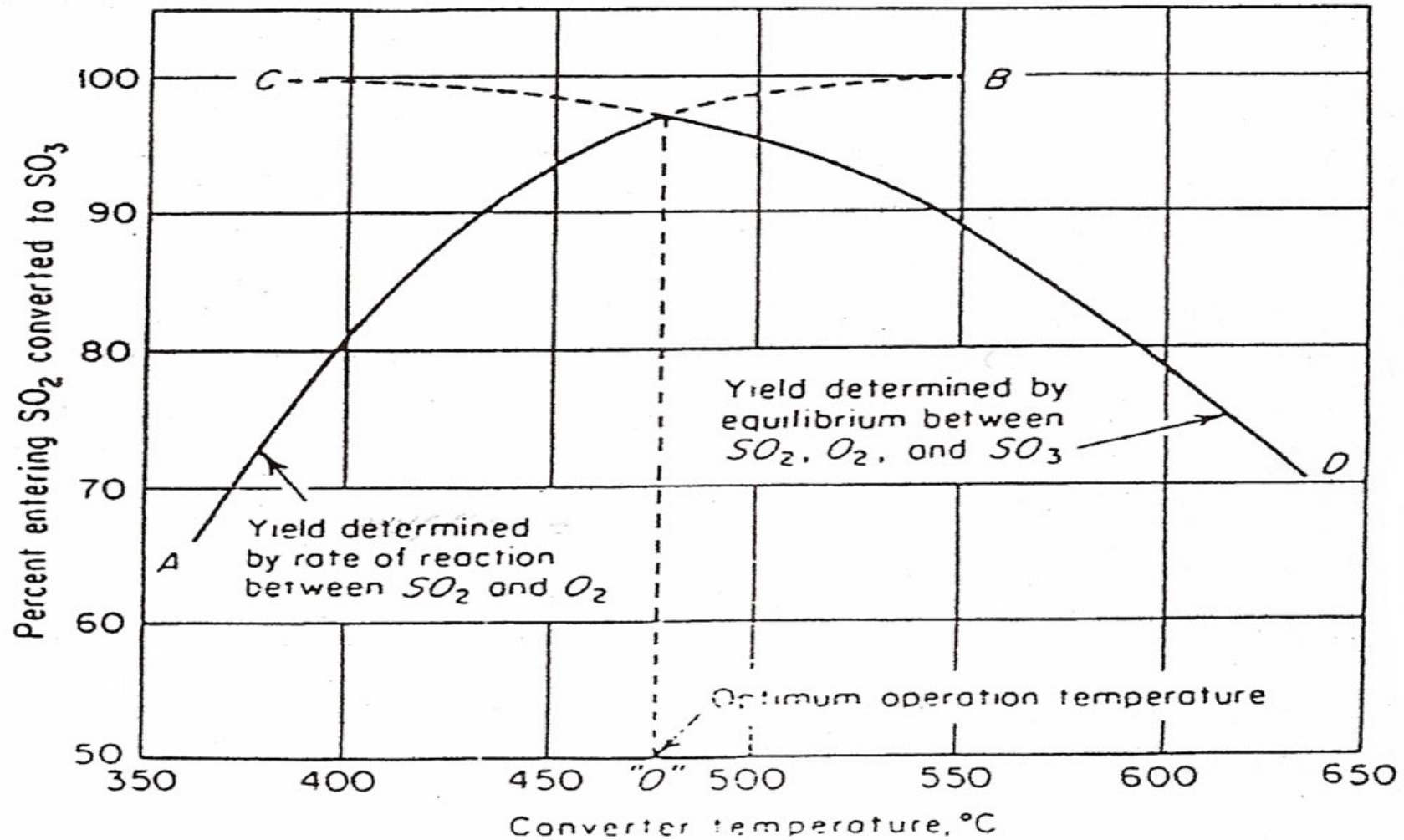

Determination of optimum operation temperature in sulfur dioxide converter.




# Introduction to Kinetics:

## Factors That Affect the Rate of Reaction

### C. Effect of Adding a Catalyst:

Catalyst does not get consumed by the reaction but it donates or absorbs electrons, hydrogen ions or hydroxide ions (see p23-4 of Experiment E4A lab manual pdf file)

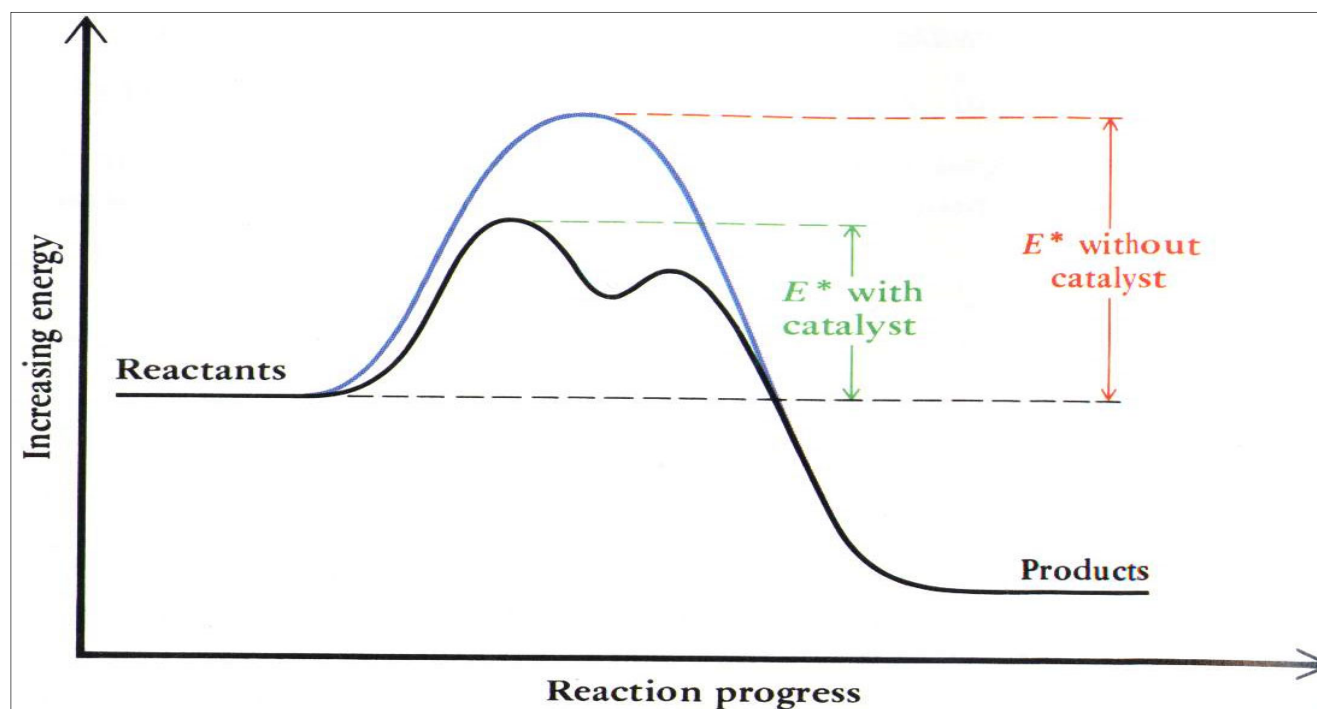

# E4 (I)

## (Two Students)

*Concentrations of reactants* →

1M HCl v.s. 6M HCl  
(dissolving eggshell)  
Dropper Plate

*Surface area* →

Fe wire v.s. Fe powder  
(w/  $\text{CuSO}_4$ )  
Two beakers @ 80 °C

*Temperature* →

Ice cold v.s. 80 °C  
Fe powder w/  $\text{CuSO}_4$

*Catalysts* →

Ice cold v.s. 80 °C  
 $\text{H}_2\text{O}_2$  decomposition w/  $\text{MnO}_2$

	A	B	C	D	E	F	G	H	I	J	K
1	<b>VC211 EXPERIMENT E4(I) DATASHEET: KINETICS FACTORS AFFECTING REACTION RATES</b>										
2	SECTION: _____			TA: _____			<b>EACH 2 STUDENTS COMPLETE ENTIRE E4(I) EXPERIMENT</b>				
3	GROUP EXPERIMENT BUT EACH TWO STUDENTS DO ENTIRE EXPERIMENT. SUBMIT INDIVIDUAL REPORT										
4	PROCEDURE PARTS →		A	A	B	B	C	C	E	E	
5			CONCENTRATION EFFECT		SURFACE AREA EFFECT		TEMPERATURE EFFECT		ADDING CATALYST EFFECT		
6	GRP	NAME	ID	EGGSHELL RXNS		COLOR CHANGE $\text{CuSO}_4$ @ $80^\circ\text{C}$		COLOR CHANGE $\text{CuSO}_4$ & Zn		$\text{H}_2\text{O}_2$ & $\text{MnO}_2$ POWDER	
7	#	Chinese		1M HCl	6M HCl	Fe Wire	Fe Powder	Ice Cold	Hot @ $80^\circ\text{C}$	$80^\circ\text{C}$ , in 30-40s	ICE AFTER 30-40s
8	NOTES			Add no more than 10 drops		Use available cylinder to measure 5 mL and the available hot water bath		Use the available hot water bath		To avoid back flash avoid direct view and place the sample small beaker inside larger beaker	
9											
10	1										
11	1										
12	1										
13	1										
14	2										
15	2										
16	2										
17	2										
18	3										
19	3										
20	3										
21	3										
22	4										

- IMPORTANT NOTES**
- Each two students will do entire E4(I)**
- Clean all glassware and rinse with distilled water**
- Use 5mL  $\text{CuSO}_4$  and hold beaker down with hand while wearing cloth gloves**
- Use hot water bath for heating to  $80^\circ\text{C}$**
- Handle hot beakers with tong or cloth glove**
- Work safely & dispose chemicals in waste container**
- Must follow chemical disposal instructions: E4(I) waste in one large beaker (no rinse water), then remove solids into its own waste container, then drain solution into inorganic waste containers, while disposing E4(II) waste into another beaker then remove stirring rod and place on top of stirrer machine pan then dispose the solution in organic waste container**

Sheet1 Sheet2 Sheet3 +

# E4 (II)

## (One Group)




$$R_a = \frac{-\Delta[I_2]}{\Delta t} = \frac{-([I_2]_{\text{final}} - [I_2]_{\text{initial}})}{\Delta t}$$

$$R_a = k[A]^m[B]^n = -\Delta([I_2]_{\text{final}} - [I_2]_{\text{init}}) / \Delta t$$



	<b>4 M acetone mL</b>	<b>H<sub>2</sub>O mL</b>	<b>1 M HCl mL</b>	<b>0.00118 M I<sub>2</sub> mL</b>	<b>Total Volume mL</b>
<b>1</b>	<b>10.0</b>	<b>20.0</b>	<b>10.0</b>	<b>10.0</b>	<b>50.0</b>
<b>2</b>	<b>x1</b>	<b>y1</b>	<b>10.0</b>	<b>z1</b>	<b>50.0</b>
<b>3</b>	<b>x2</b>	<b>y2</b>	<b>10.0</b>	<b>z2</b>	<b>50.0</b>
<b>4</b>	<b>x3 (not used above)</b>	<b>y3</b>	<b>10.0</b>	<b>z3 (not used above)</b>	<b>50.0</b>
<i><b>HCl acts as catalyst (remains unchanged at end of reaction)</b></i>					

# E4 (II)

## (One Group)

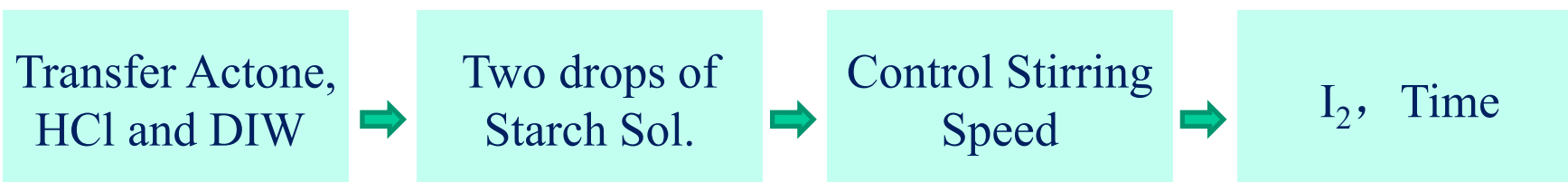

Color Change  
Enhanced

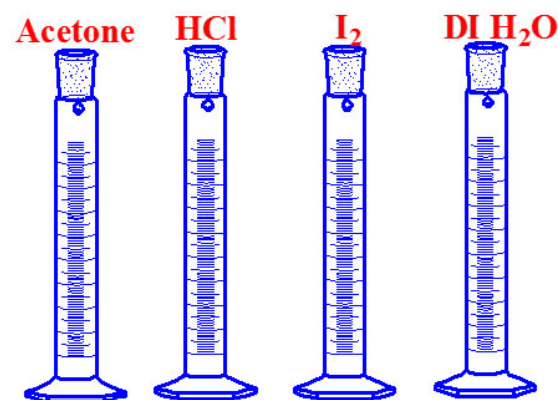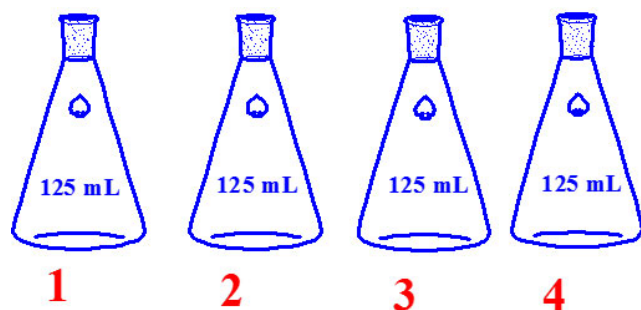


One Repeat

# VC211 EXPERIMENT E4(II) DATASHEET: DETERMINING THE RATE LAW

TA: \_\_\_\_\_

LAB ROOM: \_\_\_\_\_

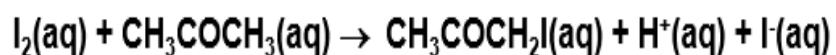

Procedure Part ----->

GROUP EXPERIMENT BUT EACH STUDENT MUST DESIGN, PREPARE & TEST ONE SAMPLE. SUBMIT INDIVIDUAL REPORTS

RATE

Compare color to clarity of blank 50 mL DI H<sub>2</sub>O

SAMPLE

1M

4M

0.00118M

INT. M

INT. M

TRI. 1

TRI. 2

AVG.

CALCULATIONS

GROUP

NAME

ID

#

HCl

Acetone

D-I Water

Iodine

Acetone

Iodine

RXN. Time

RXN. Time

RXN. Time

m n k

#

Chinese

#

mL

X mL

Y mL

Z mL

Moles/L


Moles/L

t1 (s)

t2 (s)

t<sub>avg</sub> (s)

- IMPORTANT NOTES:**
- Each student must do one reaction of entire E4(II)
- Clean all glassware and rinse with distilled water
- Use graduated cylinders to measure volumes directly
- From stock solutions, do not use beakers
- Must stir reaction without heating & splashing starting low setting
- Compare reaction color change to blank water sample – glassware over white paper
- Must not dispose stirrer rod in sink or waste container
- You will be penalized heavily if you do not follow
- Work safely & dispose chemicals in waste container
- Must follow chemical disposal instructions: E4(I) waste in one large beaker (no rinse water), then remove solids into its own waste container, then drain solution into inorganic waste container, while E4(II) waste in another beaker then remove stirring rod and place on top of stirrer machine pan then dispose the solution in organic waste container

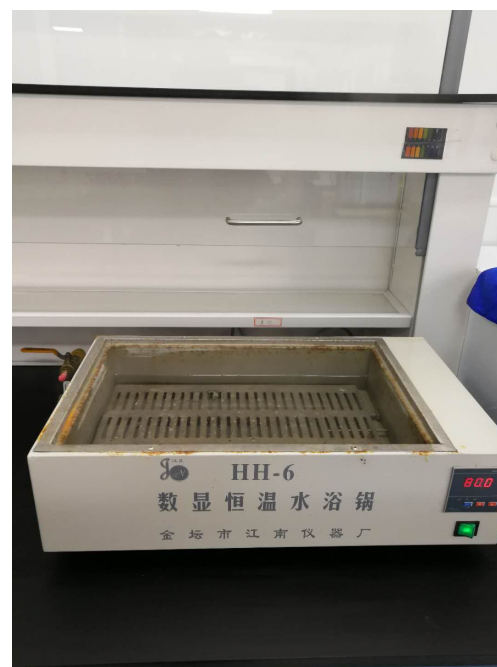

# **PART A: Introduction to Kinetics:**


## **Factors That Affect the Rate of Reaction**

**Procedure: DO NOT USE MORE THAN 10 DROPS**  
**ALWAYS HANDLE HOT BEAKERS WITH CLOTH**  
**GLOVES**

**Part A: Effect of Changing the**  
**Concentration of Reactants**

**1M HCl**

**6M HCl**

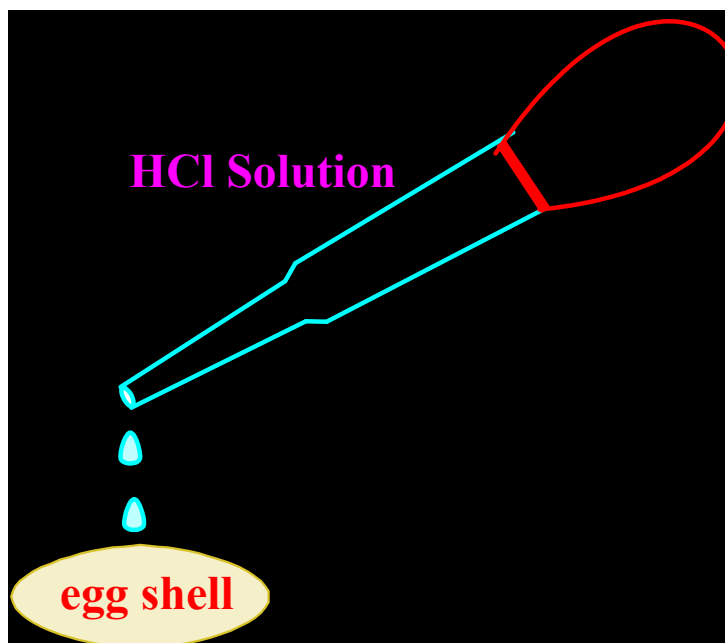




# PART B: Introduction to Kinetics: Factors That Affect the Rate of Reaction

Procedure: For heating, use the hot water bath

## Part B: Effect of Changing the Surface Area-----I


**CuSO<sub>4</sub> 5.0 ml**

IRON WIRE / IRON POWDER

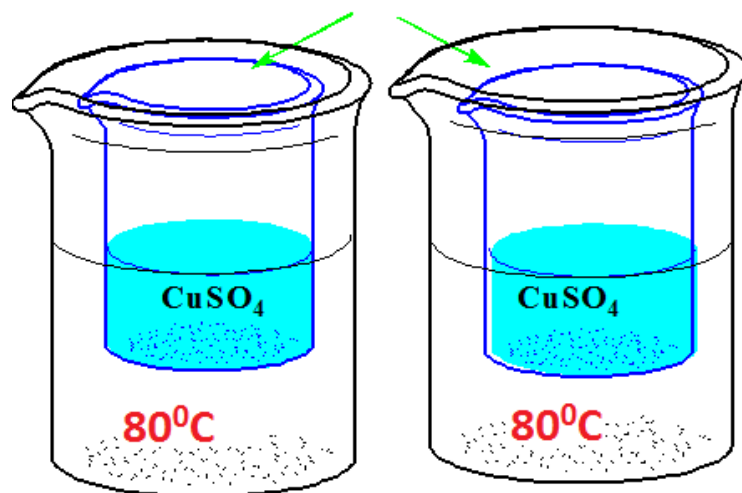


**CuSO<sub>4</sub> 5 ml**

0.2 grams  
iron wire

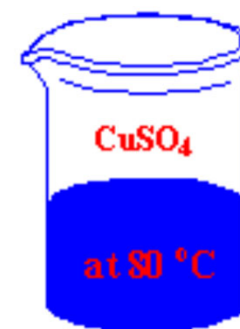

0.2 grams  
iron powder

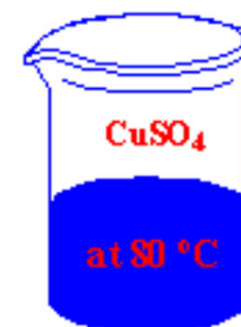

# **PART C: Introduction to Kinetics:**


## **Factors That Affect the Rate of Reaction**

**Procedure:** add 0.05 g Zn granular to  $\text{CuSO}_4$  solution

**For heating, use the hot water bath**

**Part C: Effect of Changing the Temperature-----II**

**ALWAYS HANDLE HOT BEAKERS  
WITH CLOTH GLOVES**

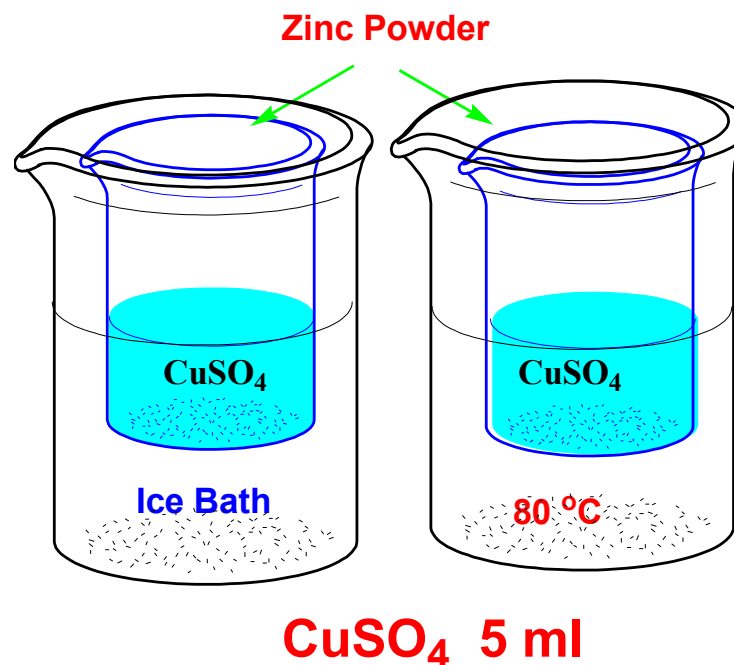

# Introduction to Kinetics:

## Factors That Affect the Rate of Reaction


### Part D: Effect of Adding a Catalyst----- Decomposition of $\text{H}_2\text{O}_2$

**Procedure: ADD 0.050g  $\text{MnO}_2$  INSTEAD OF 0.05g KI**

Reaction is splashy (& violent when more than 0.05 g solid is added).

Avoid observing directly from top when solid is added.

For heating, use the hot water bath

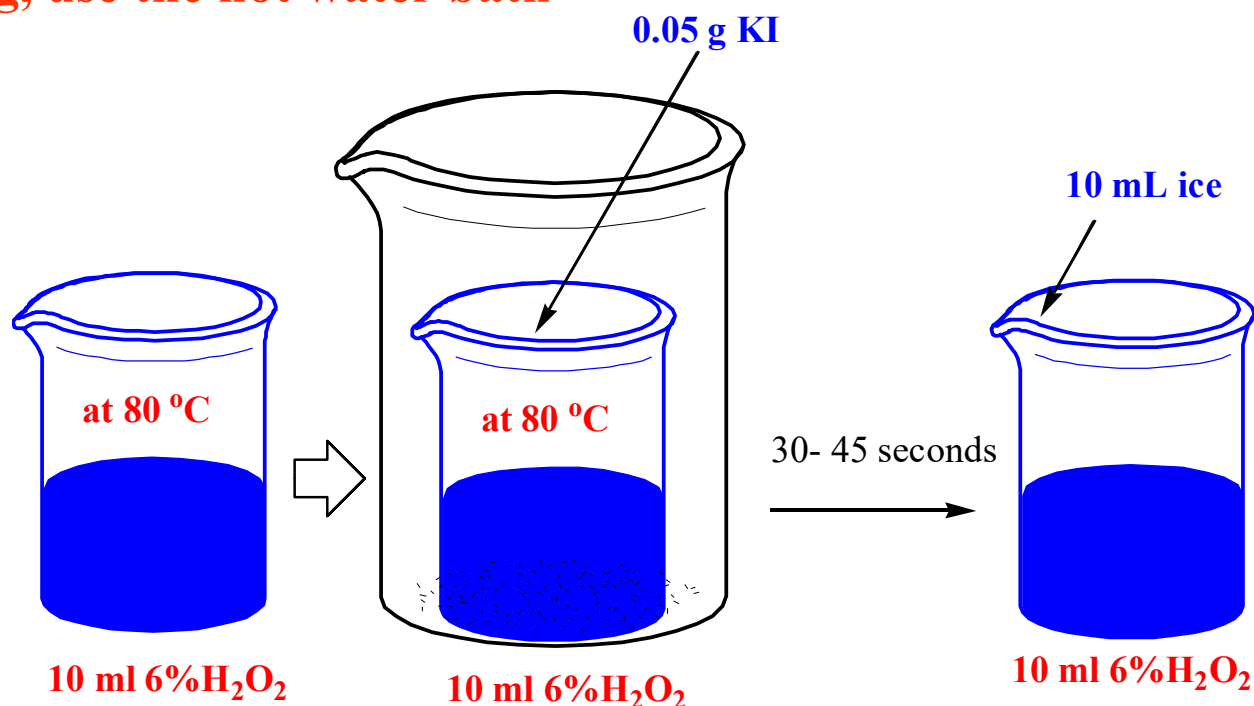


# Introduction to Kinetics:

## Factors That Affect the Rate of Reaction

**Effect of Adding a Catalyst. Catalyst does not get consumed by the reaction but it donates or absorbs electrons, hydrogen ions or hydroxide ions (see Experiment E4(I) of lab manual)**

**For heating, use the hot water bath % handle hot beakers with cloth gloves**

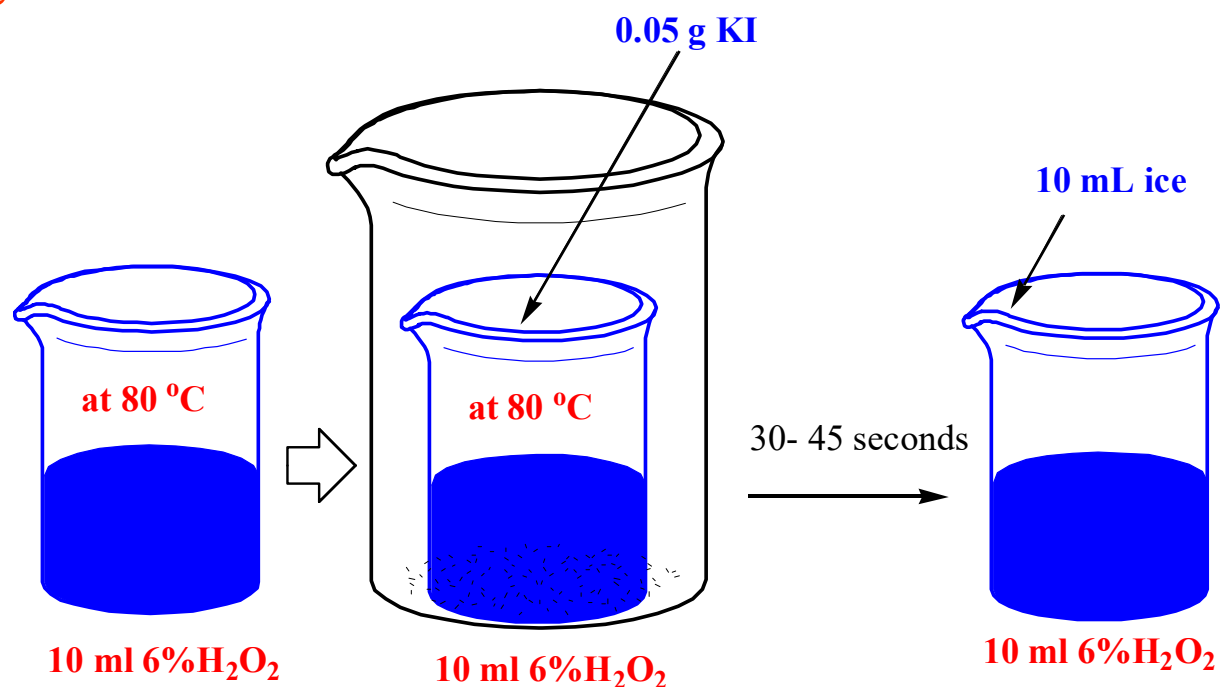

# Determining the Rate Law:

## A Kinetics Study of the Iodination of Acetone

### The Rate Law:



Reaction Rate:


$$R_A = k[A]^m[B]^n = -d[A]/dt = -\Delta[A]/\Delta t$$

**m** and **n** are determined by the experiments.

# Determining the Rate Law:

## A Kinetics Study of the Iodination of Acetone

### The Iodination of Acetone:

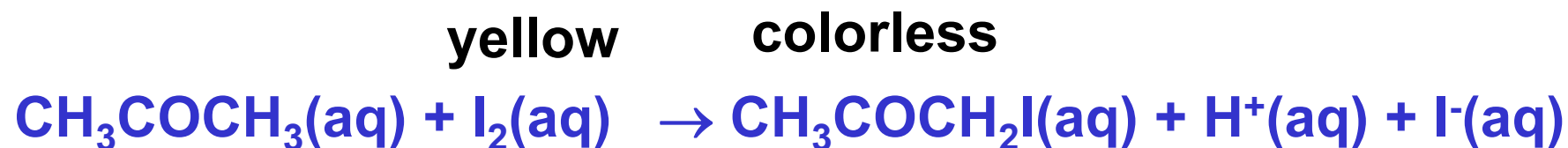

$$R_a = \frac{-\Delta[I_2]}{\Delta t} = \frac{-([I_2]_{\text{final}} - [I_2]_{\text{initial}})}{\Delta t}$$


# Determining the Rate Law:

## A Kinetics Study of the Iodination of Acetone

### Procedure:

1. Clean the glassware, use only Erlenmeyer flasks to carry the reactions & **not the beakers** (to minimize splashing effect due to stirring).

soap solution → tap water → de-ionized water

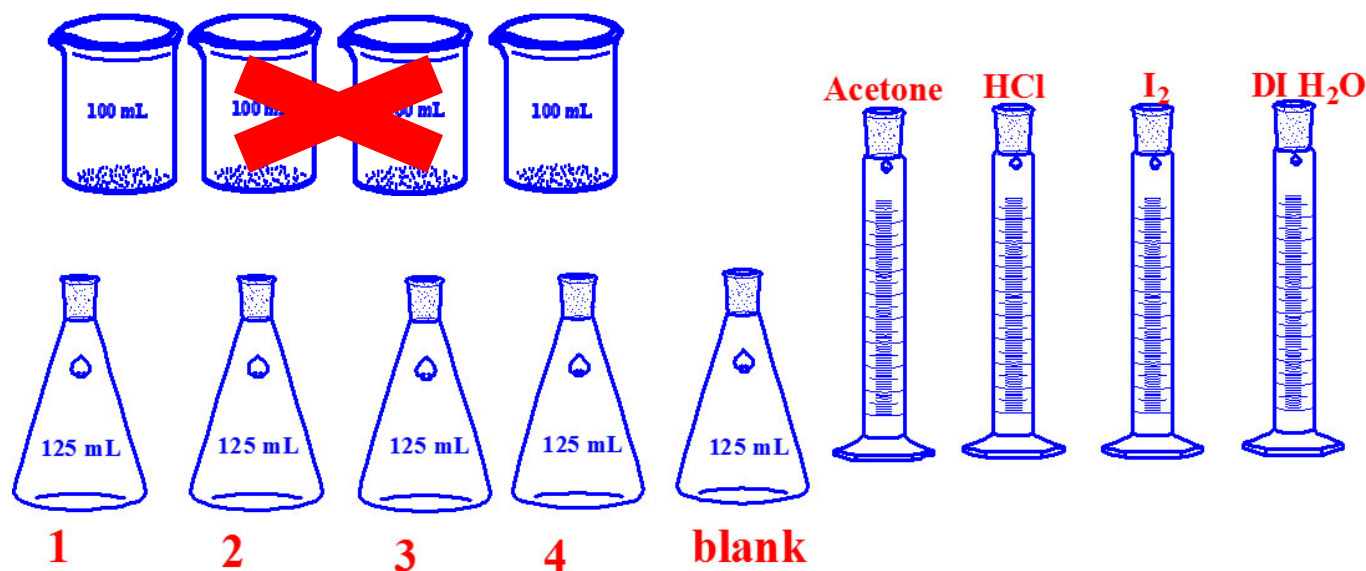


# Determining the Rate Law:

## A Kinetics Study of the Iodination of Acetone

**Procedure:** Flasks shown are 125-mL but use 250-mL or suitable Erlenmeyer flask.

2. Prepare a blank: Place over white background paper  
**50 mL of water in 250-mL Erlenmeyer flask**

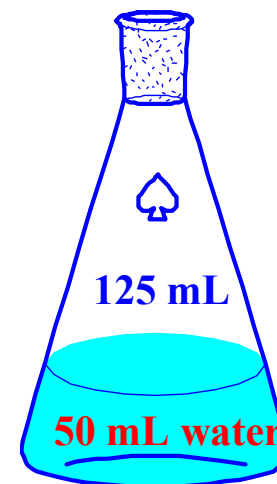




# Determining the Rate Law:

## A Kinetics Study of the Iodination of Acetone

### Procedure:

3. Prepare a mixture **without**  $I_2$ .

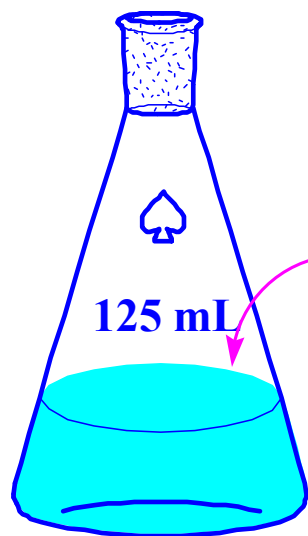

10 mL acetone  
10 mL HCl  
20 mL DI water

# Determining the Rate Law:

## A Kinetics Study of the Iodination of Acetone

**Procedure:** Place the white stirring rod inside flask

**4. Add  $I_2$  and start the timer.**

- a.  $I_2$  is added quickly with stirring**
- b. Record the time when the color of the solution changed from yellow to colorless**
- c. Record the volumes of all chemicals**

**Use white paper as background to compare REACTION color with the blank transparent water color**



# Determining the Rate Law:


## A Kinetics Study of the Iodination of Acetone

### Procedure:

5. Repeat steps 3 and 4 .

*1 time for 1 formula  
(not 2 times)*

- Calculate the percent difference between two trials.
- It is suggested to repeat until the percent difference is less than 5%, but groups do not have enough time. Therefore work diligently and. report whatever data you have.

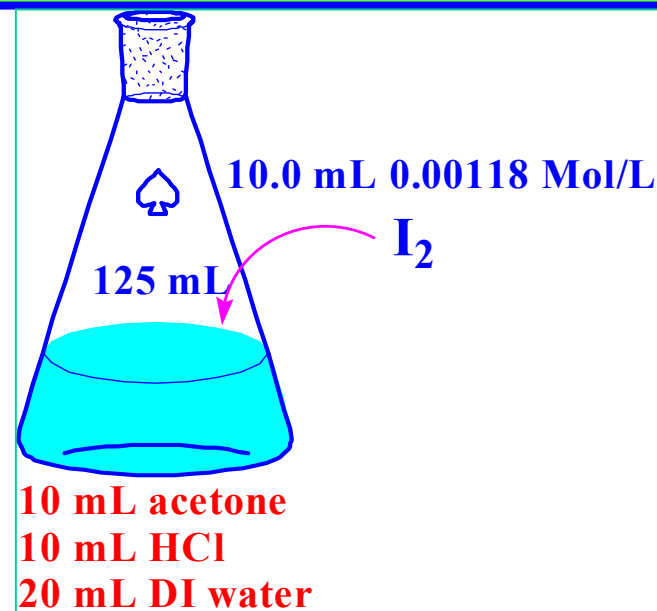


**HCl acts as a catalyst**


# Determining the Rate Law:

## A Kinetics Study of the Iodination of Acetone

### Procedure:

6. How to determine the order with respect to *Acetone* .  
*Make your proposal, and show it to your instructor.*

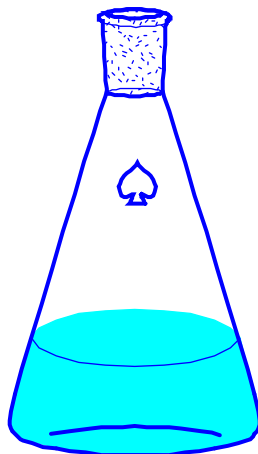

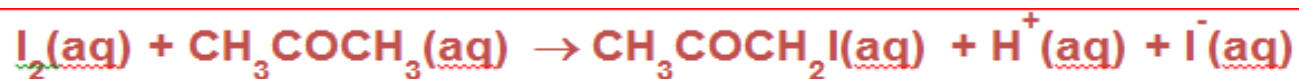

$$R_A = \frac{-\Delta[I_2]}{\Delta t} = \frac{-([I_2]_{\text{final}} - [I_2]_{\text{initial}})}{\Delta t}$$

$R_A = k[A]^m[B]^n$  Where [A] & [B] are initial conc. of I<sub>2</sub> & Acetone at -0- seconds  
 Note [I<sub>2</sub>]final = 0 when color change

Table ? : Proposed lab work to determine the reaction rate orders m & n and the reaction rate constant k.

**SEE LAST SLIDE FOR ANSWERS**

SAMPLE #	4M Acetone mL	H <sub>2</sub> O mL	1M HCl mL	0.00118M I <sub>2</sub> mL	Total Vol. mL	Initial M Acetone in 50mL Moles/L	Initial M I <sub>2</sub> in 50mL Moles/L	Trial 1 Rxn time s	Trial 2 Rxn time s	Avg. Rxn time s
1	10.0	20.0	10.0	10.0	50.0	?	?	?	?	?
2	X1	Y1	10.0	Z1	50.0	?	?	?	?	?
3	X2	Y2	10.0	Z2	50.0	?	?	?	?	?
4	X3	Y3	10.0	Z3	50.0	?	?	?	?	?

$$R_A = k[A]^m[B]^n$$

Where [A] & [B] are initial conc. of I<sub>2</sub> & acetone consecutively at time -0- seconds.

$$R_A = k[A]^m[B]^n \dots\dots\dots \text{Eq.1, but in dilution } C_{\text{conc}} V_{\text{conc}} = C_{\text{dil}} V_{\text{dil}} \dots\dots\dots \text{Eq.2}$$

$$\text{Or } C_{\text{dil}} = C_{\text{conc}} V_{\text{conc}} / V_{\text{dil}} \dots\dots\dots \text{Eq.2}$$

$$\text{But: } R_{A1} = k_1 [C_{A1}/V_{A1}]^m [C_{B1}/V_{B1}]^n \text{ or}$$

$$R_{A2} / R_{A1} = (k_2/k_1) [C_{A2}/V_{A2}]^m [C_{B2}/V_{B2}]^n / ([C_{A1}/V_{A1}]^m [C_{B1}/V_{B1}]^n) \dots\dots\dots \text{Eq.3}$$

Using Eq.2 into Eq.3 and applying the equations above at constant temperature ( $k=k_1=k_2$ ) to data of Sample #1, #2 & 3, to get the reaction order m, n & the reaction rate constant k.

Remember  $V_1 = V_2 = V_3 = V_4 = V_{\text{dil}} = 50\text{mL}$ , while  $V_{\text{conc}}$  is the proposed design volumes  $X_i$  or  $Z_i$  in the table below so the distilled water volume is  $Y_i = 50 - X_i - Z_i - 10\text{mL}$ , and the starting concentrations  $C_{\text{conc}}$  are given below for each solution.

# Method of Initial Rates

(different example than manual but same analysis)



The **method of initial rates** involves a **series of experiments** in which the **initial concentrations of some reactants are held constant** and **others are varied in convenient multiples** in order to determine the rate law for that reaction

Experiment	Initial [NO]	Initial [Cl <sub>2</sub> ]	Initial Rate, M s <sup>-1</sup>
1	0.0125 M	0.0255 M	$2.27 \times 10^{-5}$
2	0.0125 M	0.0510 M	$4.55 \times 10^{-5}$
3	0.0250 M	0.0255 M	$9.08 \times 10^{-5}$

$$\text{Rate} = k [\text{NO}]^2 [\text{Cl}_2]$$

**Using data for reaction of NO & Cl<sub>2</sub>  
we can now calculate m, n, & k**



Then:  $-\mathbf{R_a} = \mathbf{k} [\mathbf{NO}]^m [\mathbf{Cl}_2]^n$

From data in previous table:

$$2.27 \times 10^{-5} = k (0.0125)^m (0.0255)^n \dots\dots\dots \text{Eq.1}$$

$$4.55 \times 10^{-5} = k (0.0125)^m (0.0510)^n \dots\dots\dots \text{Eq.2}$$

$$9.08 \times 10^{-5} = k (0.025)^m (0.0255)^n \dots\dots\dots \text{Eq.3}$$

$$\text{Eq.2/Eq.1: } 2 = 2^n, \text{ then } n = 1$$

$$\text{Eq.3/Eq.1: } 4 = 2^m, \text{ then } m = 2$$

**The overall order of reaction = m + n = 3, first order  
wrt Cl<sub>2</sub> gas & second order wrt NO gas.**


**Using any of the equations, then k is solved:**

$$\mathbf{k} = \mathbf{5.70 \text{ M}^{-2} \text{ s}^{-1}}$$

$$-\mathbf{R_a} = \mathbf{5.70 \text{ M}^{-2} \text{ s}^{-1} [\text{NO}]^2 [\text{Cl}_2]}$$



## Example 2: REACTION ORDER

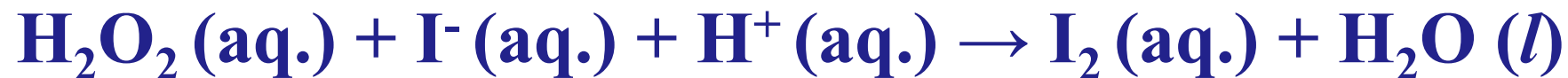

**The rate law expression for the reaction:**

**Rate of disappearance of reactants = - Rate of appearance of products.**

**I<sub>2</sub> was easily observed for color change. Then the rate of production of I<sub>2</sub> :**

$$R_{\text{I}_2} = v = k [\text{H}_2\text{O}_2]^m [\text{I}^-]^n [\text{H}^+]^p$$

**Note: Using enough experimental data we can easily determine m, n, p & k**

**DATA TABLE:** Proposed lab work to determine the reaction rate orders m & n and the reaction rate constant k.

SAMPLE #	4M Acetone mL	H <sub>2</sub> O mL	1M HCl mL	0.00118 M I <sub>2</sub> mL	Total Vol. mL	Initial M Acetone in 50mL Moles/L	Initial M I <sub>2</sub> in 50mL Moles/L	Trial 1 Rxn time s	Trial2 Rxn time s	Avg. Rxn time s
1	10.0	20.0	10.0	10.0	50.0	? 0.8	? $2.36 \times 10^{-4}$	?	?	?
2	X1 10.0	Y1 10.0	10.0	Z1 20.0	50.0	? 0.8	? $4.72 \times 10^{-4}$	?	?	?
3	X2 20.0	Y2 10.0	10.0	Z2 10.0	50.0	? 0.8	? $2.36 \times 10^{-4}$	?	?	?
4	X3 20.0	Y3 0.0	10.0	Z3 20.0	50.0	? 0.8	? $4.72 \times 10^{-4}$	?	?	?
BLANK	0	50.0	0	0	50.0	Compare sample solution color to water sample transparent color				

**Notes:**

► The total volume should be 50.0mL, so Y values should be the difference between 50.0mL and what